

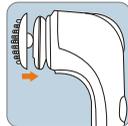
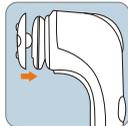


# Multi-Functional Massager CO2

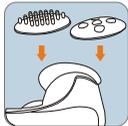

## FUNCTION&CHARACTERISTIC

- Vibration massage works together with hot pack leisure device, enjoys double from one massage.
- The streamline long handle design, it is very easy to take and operate.
- Non-step speed control, you can adjust the machine freely.
- Attaches two sets detachable massage parts.
- Intellectualized time control design, safe and comfortable.

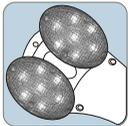

## USAGE

- Please make sure the pattern of massage stick point at 'off' before have power
- Adjust the pattern of the massage stick point at "massaging" or "massaging and hot compress function"
- Adjust the speed button of the massage stick
- If you running the machine according to the instruction, will be make your relax, don't worried about harm to your muscle.
- Only need back and forth to move, each place don't surpass 3 mins, and adjust the place and time of the Massage according to your requirement. in order to avoid damage.
- Attaches two kinds of detachable cushion (finger-pressure type and carding brushing type), install it according to the chart.
- Please depend on your enjoy, adjust the speed of the vibration.

Operation mode



Off




Massaging And  
Hot Compress Function




Speed

# WARNING

- Please do not use the product under a voltage other than that stipulated for the product.
- Please keep the controller away from water to prevent damage.
- Please do not use this product in highly humid environments such as bathrooms, swimming pools and hot springs.
- Please do not place the product in the vicinity of heaters.
- Please do not use the product if the casing or case has been damaged, as this might result in injury.
- Damaged casing or cases should be replaced as soon as possible.
- Please make sure that the power switch is off before inserting the plug into the socket.
- Please do not attempt to remove the plug with wet hands or by pulling the cable.
- Please always remember to turn the power switch to the off position and remove the plug from the socket correctly after use.
- Please do not disassemble the product if the cord, plug or any other part is damaged. Please contact your local dealer or repair station for maintenance or repair.
- Certain parts of the product may become hot following a long period of use. People that are sensitive to Heat should not use the product for long periods of time.
- The product may stop working after a long period of use to prevent internal overheating. It will work again normally once it has cooled-down.
- Please don't running the machine surpass above 20mins.
- Children under 14 years of age or people suffering from mental disorders are prohibited from using this Product without supervision.
- Pregnant and menstruating women are advised not to use this product.

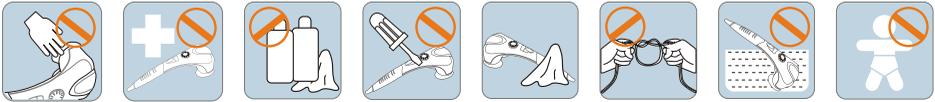

## NAME OF EACH PART

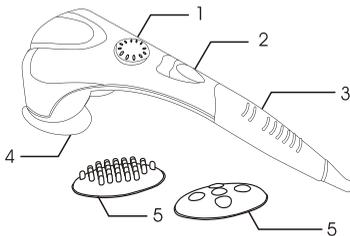

1. Adjustable speed knob
2. Power/use mode switch
3. Lengthening handle
4. Leisure hot pack device
5. Changeable massage mat

## JUGEMENT OF FAULT

Fault	Reason	Fault Remedy
Motor not working	Plug not fixed up	Fix up the plug well.
	Switch is out of position	Checking the Switch
Working too long	Thermal Protection	Rest for 5 minutes

## PRODUCT SPECIFICATIONS

Product model : C02

Name: Multi-Functional Massager

Configure:02

Input Voltage: 220-240V~ 50Hz/60Hz  110-120V~ 60Hz

Rated power: 28W

Safety structure: